


CT Tibial Plateau

Application Examples: fracture	
Scan Protocol (See for Technical Factors)	Knee

Reformations: Coronal and sagittal MPRs. 3D VRs upon request.

1. *2x2 mm coronal and sagittal MPRs* orthogonal to affected knee. Make axial MPR data set if needed.

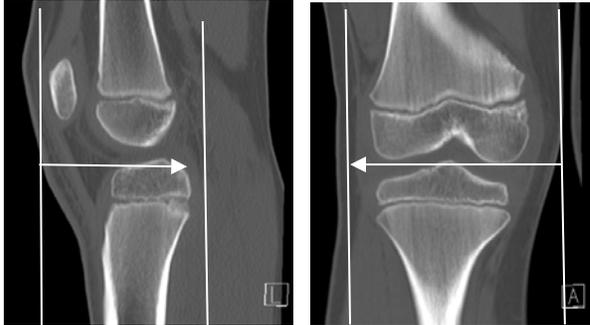

2. *3D VR of entire knee.* 360 degree rotation and image number set at 36.

Series Description: **3D**
Series Number: **1001**



3. *3D VR of tibial surface.* Start with tibia in AP position and rotate 360 degrees then tilt down and capture images every 10 degrees, until just past tibial surface.

Series Description: **3D**
Series Number: **1002**

