


## Soft Tissue Neck

Application Examples: neck cancer

Sequence	FoV	Thickness	Gap
COR T1	25CM	4MM	1MM
COR T2 IDEAL	25CM	4MM	1MM
AX T2 IDEAL	22CM	4MM	1MM
AX T1	22CM	4MM	1MM
AX DWI SOFT TISSUE NECK ( <i>VALUES B-0, B-800</i> )	22CM	4MM	1MM
INJECT CONTRAST (1mmol/Kg)			
COR T1 IDEAL WITH GAD	25CM	4MM	1MM
AX T1 IDEAL WITH GAD	22CM	4MM	1MM
SAG T1 IDEAL WITH GAD	25CM	4MM	1MM

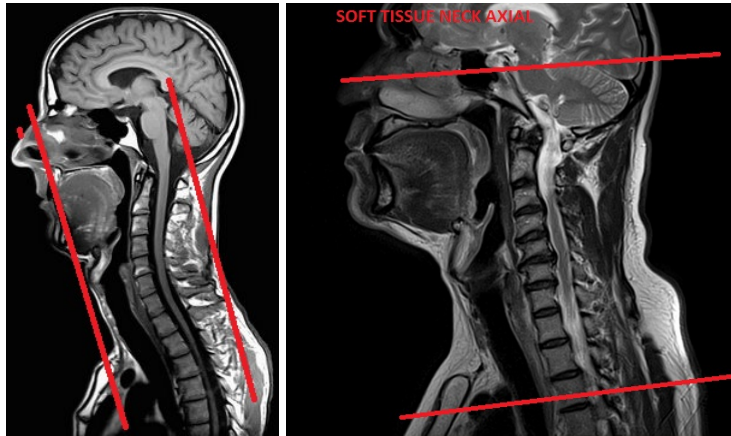

COR Coverage

AXIAL Coverage